




A Professional POCT Source Factory

Committed to Delivering Accurate, Convenient and Stable Products.

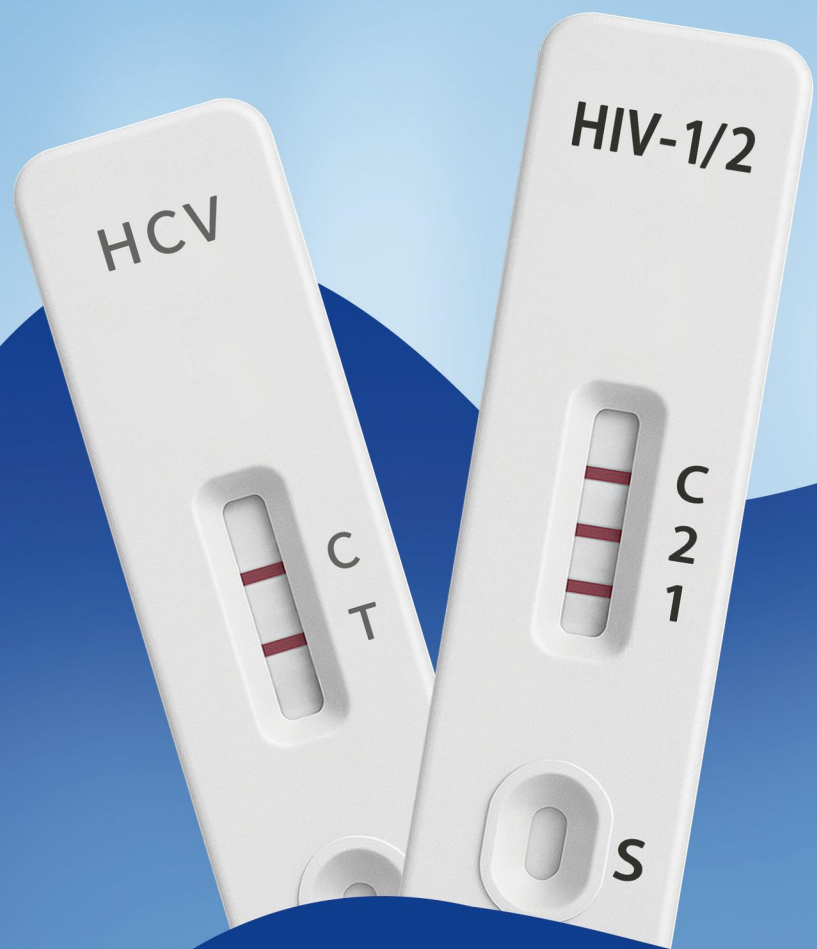

Product Catalogue

IVD MANUFACTURER SINCE 2003

NewScen Coast Bio-Pharmaceutical Co., Ltd.

Contact NewScen Team :

Add : No. 65, 6th Street, Tianjin TEDA, 300457, Tianjin, P.R.China

Web : www.newscenestestkit.com (Uncut Sheet)

Web : www.newscenbiotech.com (Rapid Test)

Email : export@newscen.com



Uncut Sheet



Rapid Test



About NewScen

We are a leading manufacturer of uncut sheet rapid tests for a range of applications, including infectious diseases, tropical (mosquito transmitted) disease, drugs of abuse, tumor marker, hormone test, and more. With over 20 years of experience in the IVD industry, we have gained a reputation for accuracy, reliability, and ease of use.

At NewScen, we are committed to providing high-quality products that meet the evolving needs of the healthcare industry. Our research and development team is dedicated to developing new and improved products, and our manufacturing facilities are equipped with state-of-the-art equipment and staffed by experienced professionals.



R & D Center




Production




Workshop




QC Center

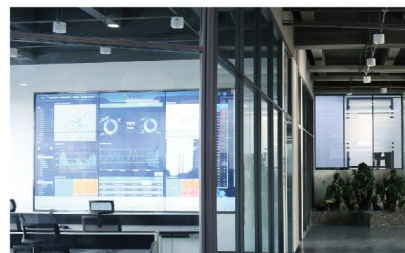


Remote Center

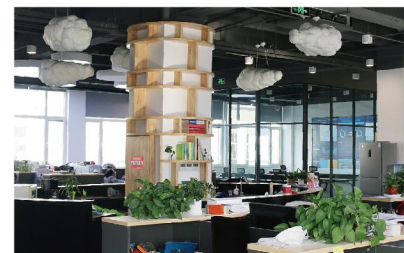

Office Center

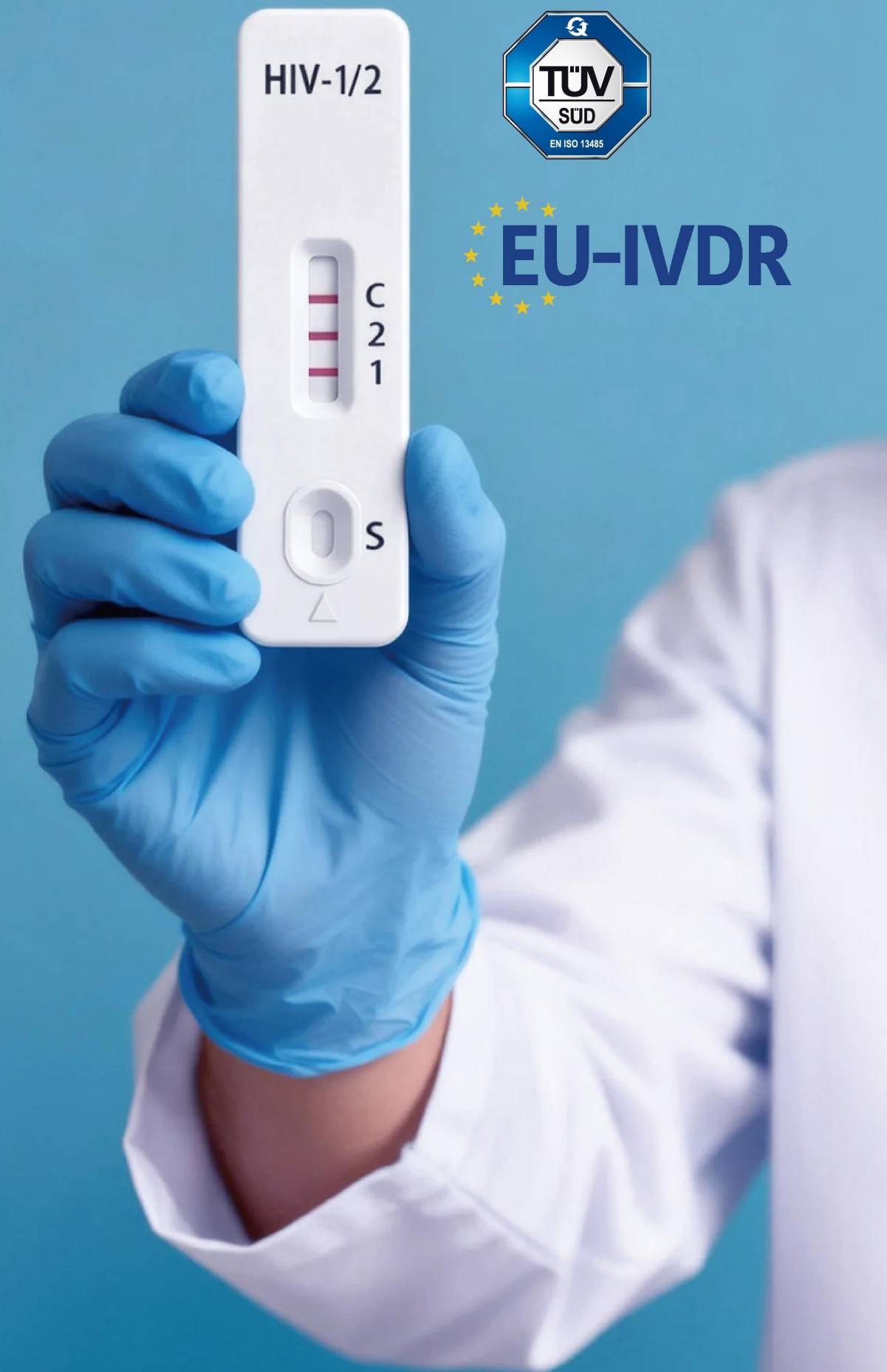

NewScen Company

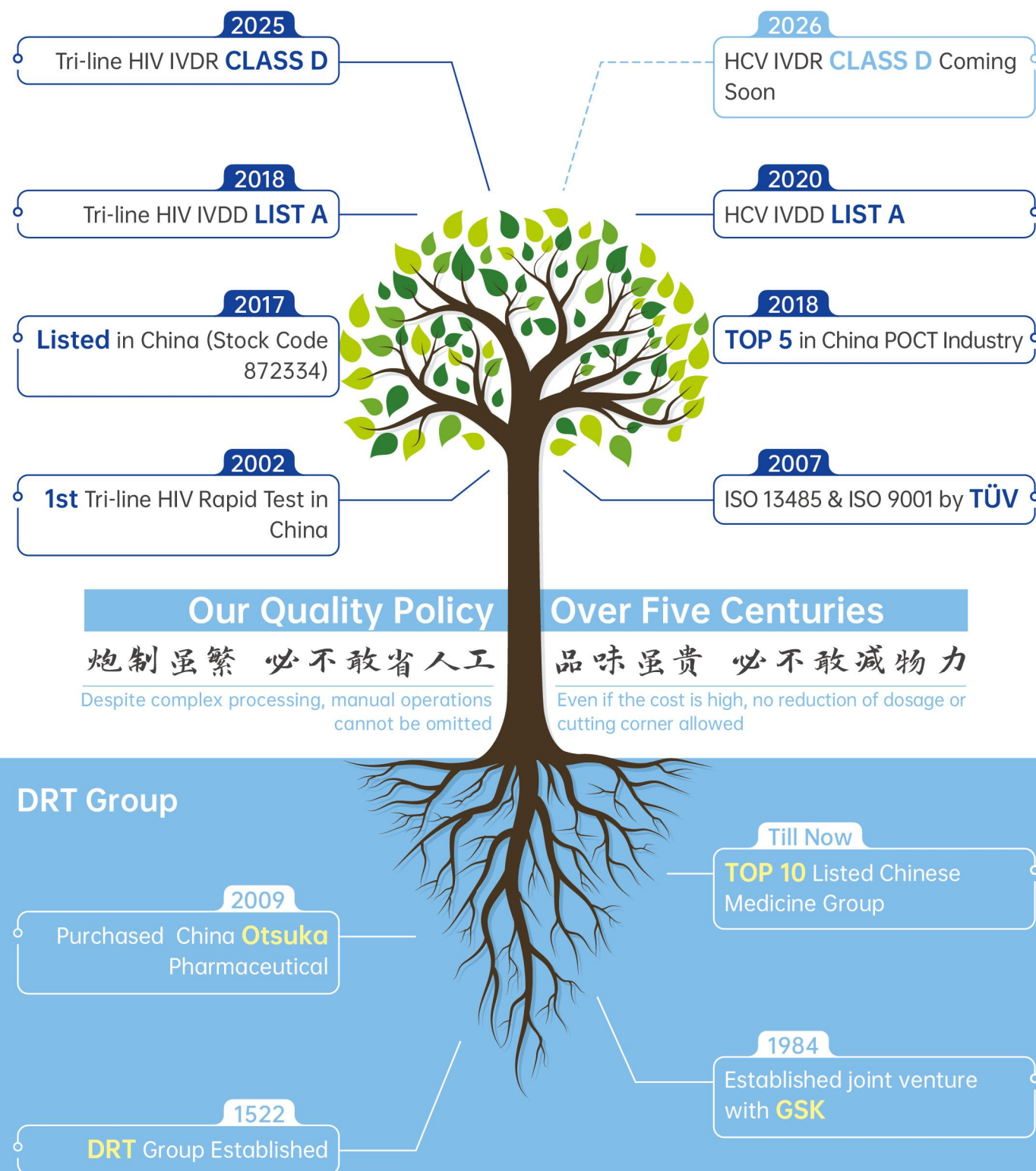

Manufacturer:

NewScen Coast Bio-Pharmaceutical Co., Ltd.

Management System:

ISO 9001 : 2015 , ISO 13485 : 2016

Detection Type Method:

Colloidal Gold

Categories:

STD, Hormone, Iron Deficiency, Tropical Disease, DOA, Tumor Marker, Respiratory Disease, Metabolism, Cardiac Marker, Inflammatory Marker, Renal Injury

Formats:

Cassette, Strip, Midstream, Uncut Sheet For Rapid Diagnostic Test Kits

Versions:

Professional, Self-Testing, OEM Options

Market:

100+ Countries and Regions Globally

Factory Tour




Display of Workshop




Components Preparations




Uncut Sheet Assembling

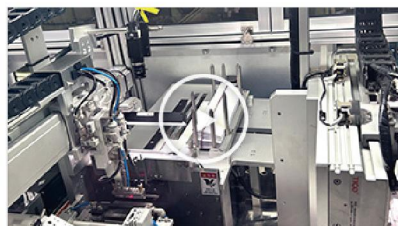

Uncut Sheet Cutting

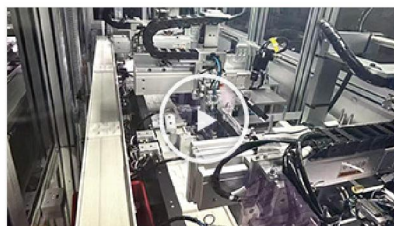

Cassette Format Processing

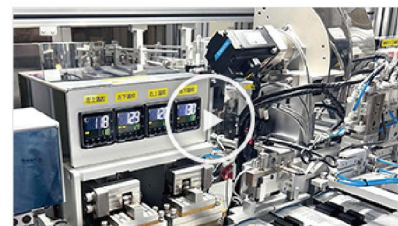

Foil Pouch Sealing



External Packing



QC Inspection of Finished Product



Warehousing for Shipping

Please scan the QR code here for NewScen factory tour.



For our colloidal gold workshop

Production Capacity: 500,000 pcs / Day



NewScen Certificates

105 Independent Intellectual Property Rights

69 NMPA Medical Device Certificates
Class III * 11 Class II * 54 Class I * 4

42 CE Certificates
CLASS D * 1 LIST A * 1 Self-Testing * 5 Others * 35

Certifications by International Authorities



Authorized Invention Patent Certificates



Colloidal Gold Test Kit

Reliable Solutions for Accurate Diagnostics







NEWSCEN HOT SALE PRODUCT LIST (2026 VERSION)

Product Name	Format					Specimen				Application		Approval
	Cassette	Strip	Midstream	Bottle	Uncut Sheet	Serum/Plasma /Whole Blood	Urine	Stool	Throat Swab/ Nasal Swab	Professional	Self Testing	
Infectious Disease Test Kit												
HIV 1/2 3 lines	√	√		√	√	√				√		NMPA / CE by TÜV / MHRA
HCV	√	√		√	√	√				√		NMPA / CE by TÜV / MHRA
HBsAg	√	√		√	√	√						NMPA
Syphilis (TP)	√	√		√	√	√						NMPA / CE others
H. Pylori (Ab)	√	√		√	√	√						NMPA / CE others
H. Pylori (Ag)	√	√		√	√			√				NMPA / CE others
TB	√			√	√	S/P						NMPA / CE others
Dengue IgG & IgM	√			√	√	√						CE others
Dengue NS1 Antigen	√			√	√	√						CE others
Malaria Pf/Pan	√			√	√	WB						CE others
Zika Virus IgG/IgM	√			√	√	√						CE others
Zika NS1 Antigen	√			√	√	√						CE others
Chikungunya IgG/IgM	√			√	√	√						CE others
Iron Deficiency Test Kit												
SF	√			√	√	√					√	NMPA / CE by TÜV / MHRA
Metabolism Test Kit												
25-OH-VD	√			√	√	√						NMPA
Hormone Test Kit												
HCG	√	√	√	√	√		√				√	NMPA / CE by TÜV / MHRA
Ovulation (LH)		√	√	√	√		√				√	NMPA / CE by TÜV / MHRA
Menopause (FSH)	√	√	√	√	√		√				√	NMPA / CE by TÜV / MHRA
DOA Test Kit												
Morphine (MOR)	√	√		√	√		√					NMPA
Methamphetamine (MET)	√	√		√	√		√					NMPA
5-IN-1 Multiple Drug Test		√			√		√					CE others
Tumor Marker Test Kit												
FOB	√	√		√	√			√			√	NMPA / CE by TÜV / MHRA
Transferrin (TF)	√			√	√			√				NMPA / CE others
AFP	√	√		√	√	√						NMPA / MHRA
CEA	√			√	√	√						MHRA
Others												
Reagent Strips for Urinalysis (chemistry method)		√		√			√					NMPA
Coronavirus Test Kit												
COVID-19 Antigen	√				√				√			CE others
COVID-19 Ag. & Influenza Combo	√				√				√			CE others

★ Semi Quantify the Test Results with NEWSCEN Solution, More Precise, More Flexible and More Cost Efficient.

Product Name	Format			Specimen		Application		Approval
	Cassette	Strip	Uncut Sheet	Serum/Plasma /Whole Blood	Urine	Semi-Quantitative	Quantitative	
Cardiac Marker Test Kit								
cTnI	√		√	√		√	√	NMPA / CE others
D-Dimer	√		√	P/WB		√	√	NMPA / CE others
H-FABP	√		√	√		√	√	NMPA
NT-proBNP	√		√	√		√	√	NMPA
Inflammatory Test Kit								
CRP	√		√	√		√	√	NMPA / CE others
PCT			√	√		√	√	NMPA / CE others
Renal Injury Test Kit								
mAlb	√	√	√		√	√	√	NMPA / CE others

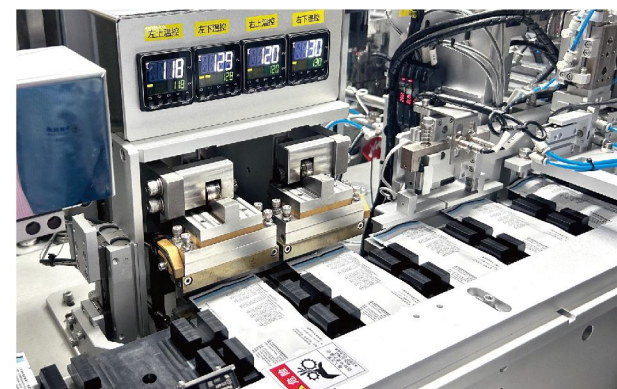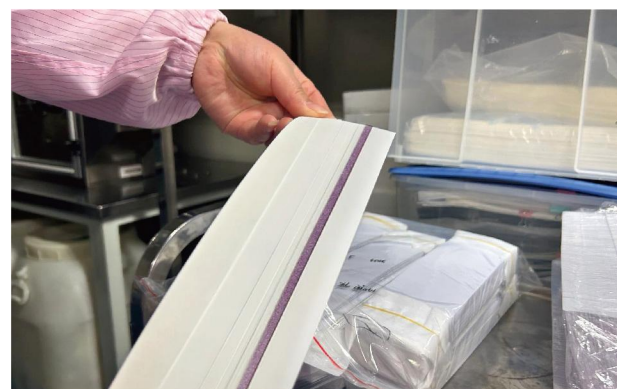



HIV (1+2) Antibody Rapid Test

Performance Advantage



IVDR Clinical Evaluation Result

Sen : **100%** ✓
Spe : **99.82%** ✓

WHO Performance Evaluation Result

Sen : **100%** ✓
Spe : **100%** ✓


Certification Threshold

TÜV SÜD high-end certification,
the highest quality endorsement



Technical Superiority

tri-lines, distinguish infection types

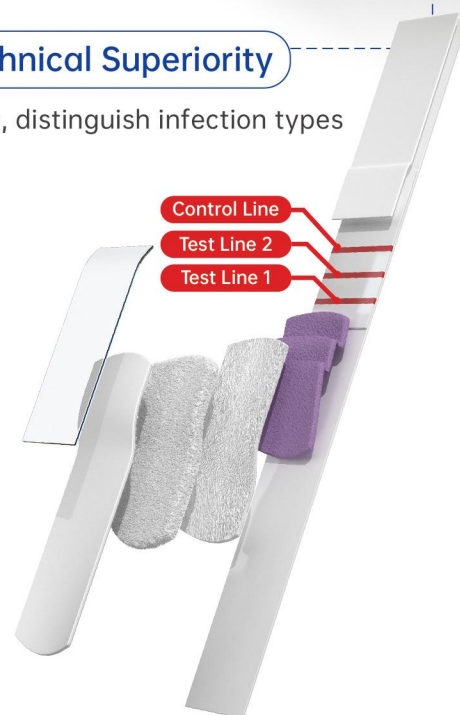

Window Period

7 days,
no competitor early than 7 days

Diversified Specification

single pack test kit available too

Market Advantage

won multiple CDC million level bids
in China

Registration Support

professional registration team,
accelerate IVDR certification cycle

Test Procedure In 3 Simple Steps



Kit Components

Components	Quantity (pcs)				
	1 Test (HIV _s -29200111000)	2 Tests (HIV _s -29200201000)	5 Tests (HIV _s -29200501000)	25 Tests (HIV _s -29202511000)	40 Tests (HIV _s -29204017000)
Test cassette	1	2	5	25	40
Desiccant	1	2	5	25	40
Dropper	1	2	5	25	40
Diluent buffer	1×1mL	1×1mL	1×1mL	2×5mL	2×5mL
Sterile lancet	1	2	5	Not provided	Not provided
Alcohol pad	1	2	5	Not provided	Not provided
Instruction for use	1	1	1	1	1

Clinical Performance

		Specimens		Sensitivity
Population	Specificity	HIV-1 Ab positive (serum/plasma)		100%
Blood donor	100%	HIV-1 Ab positive (whole blood/plasma couples)		100%
Whole blood donation	100%	non-B-subtypes HIV-1 Ab positive (plasma)		100%
Hospitalized patient	98.5%	HIV-2 Ab positive (serum/plasma)		100%
Total	99.82%	Total		100%



HCV Antibody Rapid Test

Performance Advantage



IVDR Clinical Evaluation Result

Sen : **99.76%** ✓
Spe : **100%** ✓

WHO Performance Evaluation Result

Sen : **100%** ✓
Spe : **98.7%** ✓

Certification Threshold

TÜV SÜD high-end certification,
the highest quality endorsement
(IVDR Coming Soon)




EU-IVDR



Technical Superiority

simultaneously detect
IgG, IgM, IgA

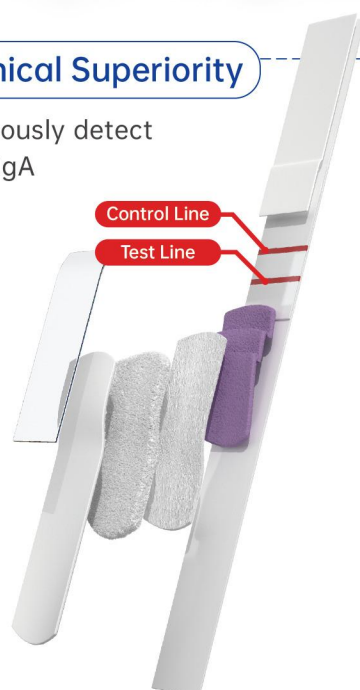

Window Period

7 days,
earlier than the benchmark assay

Diversified Specification

single pack test kit available too


Market Advantage

won million level bid in EU and
multiple CDC million level bids in China

Registration Support

professional registration team,
accelerate IVDR certification cycle

Test Procedure In 3 Simple Steps

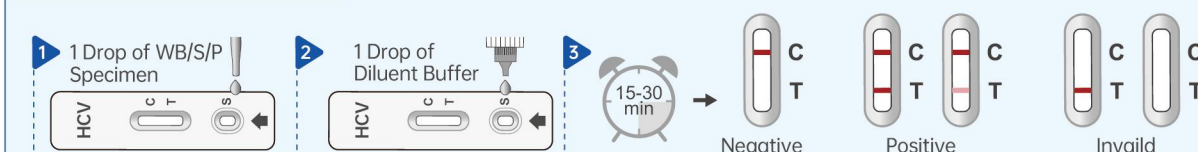

Kit Components

Components	Quantity (pcs)				
	40 Tests (HCV-29204013000)	40 Tests (HCV-29204011000)	25 Tests (HCV-29202501000)	25 Tests (HCV-29202502000)	1 Test (HCV-29200103000)
Test cassette	40	40	25	25	1
Desiccant	40	40	25	25	1
Dropper	40	40	25	25	1
Diluent buffer	2×5mL	2×5mL	2×5mL	2×5mL	1×1mL
Sterile lancet	40	Not provided	25	Not provided	1
Alcohol pad	40	Not provided	25	Not provided	1
Instruction for use	1	1	1	1	1

Clinical Performance

Population	Specificity	Specimens	Sensitivity
Blood donor (plasma)	100%	Anti-HCV positive serum	99.54%
Whole blood donation	100%	Anti-HCV positive whole blood/ plasma couples	100%
Hospitalized patient specimens	100%	Anti-HCV positive samples with known genotype	100%
Total	100%	Total	99.76%

Serum Ferritin Test Kit (SF)

Check Iron Deficiency in Minutes
Simple, Fast, Reliable!



Product Advantages

- *Identify anaemia symptoms earlier
- *Rapid 15–20minute interpretation of results
- *Accurate and stable colloidal gold method
- *High sensitivity with clear color development
- *Home-use friendly – hospital-free, time-saving, and cost-effective



What are the effects of low ferritin?



Dizziness



Palpitations



Fatigue & weakness



Dry skin



Decreased immunity



Pale complexion



Loss of appetite



Rapid heartbeat

Components



1 Test / Kit



1 × Foil Pouch with Test Cassette and Desiccant



1 × Disposable Quantitative Sampling Dropper



1 × Diluent Buffer



1 × Sterile Lancet



1 × Alcohol Pad



1 × Color Card



1 × Instruction for Use



1 × Quick Manual (Optional)

Product	Ref. No.	Pack Size
Serum Ferritin Test Kit (SF)	SF-B-29304004000	40 Tests / Kit
	SF-B-29300101000	1 Test / Kit

Storage : 4°C - 30°C · Shelf Life : 24 Months

Procedures

I. SAMPLE COLLECTION



Disinfect




Puncture




Absorb

II. TEST PROCEDURE

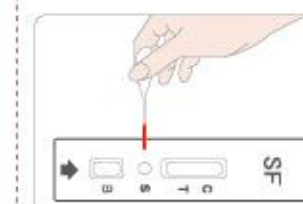

1 Drop of Blood

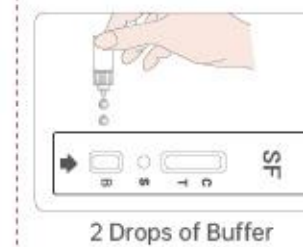


2 Drops of Buffer

III. RESULT INTERPRETATION



Within
15 - 20
minutes

Compare the Test Line Color Intensity with the Color Card

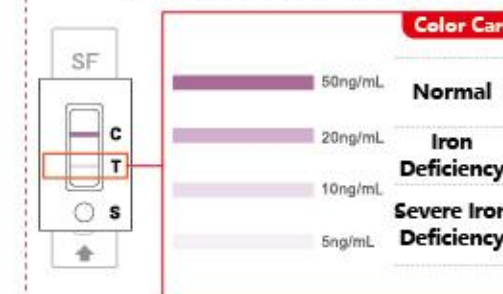

One-Step Fecal Occult Blood (FOB) Diagnostic Kit



Sensitivity: **95.7%** ✓
Specificity: **99.2%** ✓
Total Accuracy: **97.5%** ✓


Product Advantages

CE-marked self-testing kit for early detection of gastrointestinal (GI) bleeding
Combines immune and chemical methods for enhanced accuracy and wider detection range
Detects both upper and lower GI tract bleeding with high sensitivity and specificity
Patented one-step design for easy and reliable combo test results
Rapid results in 1–5 minutes, with no cross-reaction or Hook Effect interference
No need for extra equipment; non-invasive, hygienic, and user-friendly
Stable dry reagent, room temperature storage, perfect for **at-home use**

Common Symptom




Assay Procedure

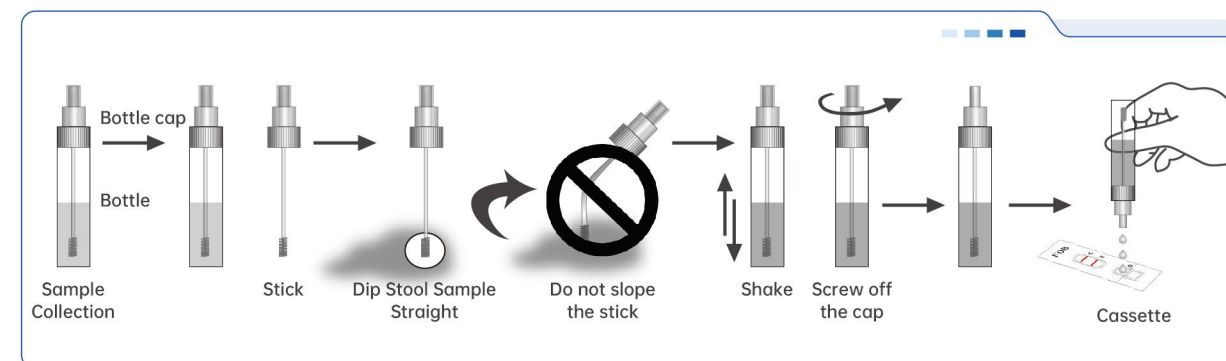

Result Interpretation

Color card (Chemical method) Observe within 1 min					
Immune Method Test Result	Positive (+)	Negative (-)	Positive (+)	Negative (-)	Negative (-)
Chemical Method Test Result	Positive (+) 2µg/mL 20µg/mL 200µg/mL 2000µg/mL	Negative (-) 0µg/mL	Negative (-) 0µg/mL	Positive (+) ≤ 200 µg/ml 20µg/mL	Positive (+) ≥ 200 µg/ml 2000µg/mL
Final Interpretation	Positive (+)	Negative (-)	Positive (+)	Negative (-)	Sample diluted 100 times & retest

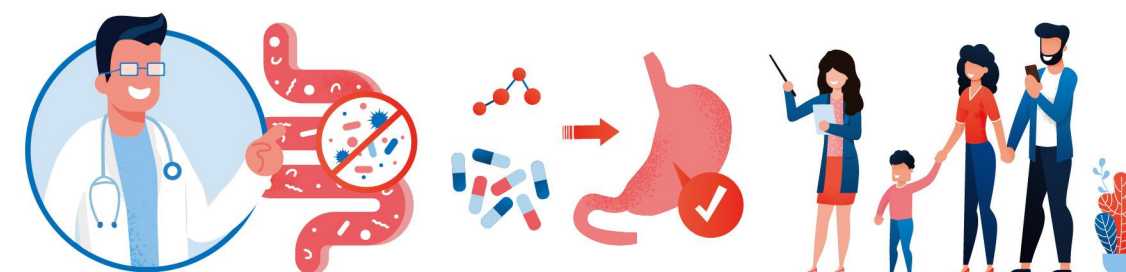

Human Chorionic Gonadotropin (HCG) tests



Product Advantages

- Clear and easy-to-read results – early pregnancy visible at a glance
- Detects pregnancy as early as 3–7 days after conception, including potential ectopic pregnancy
- High accuracy and sensitivity, detecting HCG levels as low as 25 mIU/mL
- Simple home self-test with fast results in just minutes
- Individually sealed, clean and hygienic packaging for safe and convenient use

Luteinizing Hormone (LH) tests (Ovulation tests)



Product Advantages

- Cut-off value: 20 mIU/mL
- Accurate Timing: Detects the LH surge within hours of ovulation
- Test Principle: Peak-to-fade pattern confirms true ovulation
- Easy to Use: Simple one-step home self-test
- Smart Detection: Helps identify low peak ovulation and lock in the best time to conceive

Follicle Stimulating Hormone (FSH) tests (Menopause tests)



Product Advantages

- Suitable for women of different ages to assess ovarian function and fertility status
- Fast and accurate results using colloidal gold technology
- Convenient home self-test – private, discreet, and user-friendly
- Can be tested anytime, regardless of menstrual cycle
- High sensitivity (as low as 25 mIU/mL), hospital-free and cost-effective



H. Pylori (HP) Antigen Rapid test



Product Advantages

- Non-invasive detection of H. Pylori infection through stool samples
- Rapid results in 10–15 minutes for quick clinical decisions
- High accuracy using colloidal gold immunoassay technology
- Easy to operate with clear visual interpretation
- Suitable for home or clinical use – no need for endoscopy or biopsy

H. Pylori (HP) Antibody Rapid test



Product Advantages

- Detects H. Pylori infection accurately through antibody screening
- Quick results within 10–15 minutes for early intervention
- Utilizes colloidal gold method – stable, reliable, and easy to interpret
- Non-invasive testing using a fingerstick whole blood sample
- Ideal for home or clinic use – simple, fast, and no lab required

Significance of Combined Detection of HPAg and HPAb

Sample

Non-invasive, multi sample types for more reliable early diagnosis of H.Pylori. Convenient collection of stool specimen for HPAg and easy sampling of finger tip whole blood specimen for HPAb (S/P/WB sample applicable for professional use).

Operation

Easy operation for professional use in labs, hospitals and community health service centers (applicable for S/P/WB and stool sample) as well as the self sampling users at home for prognostic monitoring regularly while protecting personal privacy.

Interpretation

Combo detection of HPAg and HPAb offers 4 possible combinations of test results for prompt and clear interpretation of the H.Pylori infection and/or recovery status simultaneously within 10-15 minutes and provides reliable reference to timely and proper treatment.



Treponema Pallidum (TP) Antibody Rapid Test Kit



Product Advantages

- Cut-off value: 1NCU/mL. High sensitivity, specificity, accuracy and stability for clinical testing
- Patented technology, double antigen sandwich method for rapid screening of syphilis
- Reliable test result in consistent with that of the TPPA method detection kit
- Applicable for fingerstick whole blood, serum and plasma sample
- One-step method, easy operation, no equipments required
- Suitable for POCT testing, flexible application scenarios

Tropical Disease Rapid Test

MAR: Malaria (Pf/Pan)

Product Advantages

- Comprehensive and rapid detection of all 4 species (P.f, P.v, P.m, P.o)
- Differential diagnosis between P.f and other Plasmodium species
- Immediate and reliable double confirmation for P.f infection
- Room temperature storage (2-30°C), long shelf life (24 months)
- Applicable for both capillary and venous whole blood sample
- Easy and simple operation, suitable for point of care testing



DEN: Dengue (IgG/IgM), Dengue (NS1 Ag)

Product Advantages

IgG/IgM:

- Clear and fast differential diagnosis between primary and secondary infection
- Micro sampling (5μL), applicable for serum/plasma/whole blood sample

NS1 Ag:

- Double antibody sandwich method for early diagnosis of dengue infection
- Micro sampling (10μL), applicable for serum/plasma/whole blood sample

IgG/IgM & NS1 Ag:

- Easy operation, ideal for point of care testing in resource limited settings
- Test result can be shown within 1 minute for strong positive sample
- Room temperature storage (2-30°C), shelf life up to 24 months



ZIKA: Zika (IgG/IgM), Zika (NS1 Ag)

Product Advantages

- Differential diagnosis between primary and secondary infection
- Applicable for serum/plasma/whole blood sample
- High accuracy, clear visual interpretation
- Ideal for resource limited settings
- Room temperature storage



CHIK: Chikungunya (IgG/IgM)

Product Advantages

- Differential diagnosis between primary and secondary infection
- Applicable for serum/plasma/whole blood sample
- Convenient operation, suitable for POCT
- Room temperature storage



RDTs Series

Drugs of Abuse Rapid Test

Customization

Customized OEM service for DOA rapid test in multiple formats (Multi Panel/Cassette/Strip/Cup/Uncut Sheet)

Parameters

20+ parameters for desirable combinations of multi panel test (Urine/Saliva Sample, Regular/Customized Sensitivity)

One Step Multiple Drug Test Kit (MOR/MET/KET/THC/MDMA)



Lateral Flow Assay UNCUT SHEET

Custom Samples Tailored Solutions


Beyond Typical OEM Cooperation

Start from Uncut Sheet

The uncut sheet is consisted of **key components** as below :

- Sample Pad
- Colloidal Gold Conjugate Pad
- NC Membrane
- Absorbent Pad
- PVC Backing Sheet

Upon your receiving of the pre-assembled uncut sheet, it is strongly recommended a 3rd QC inspection procedure be properly arranged **prior to** further operation or processing.

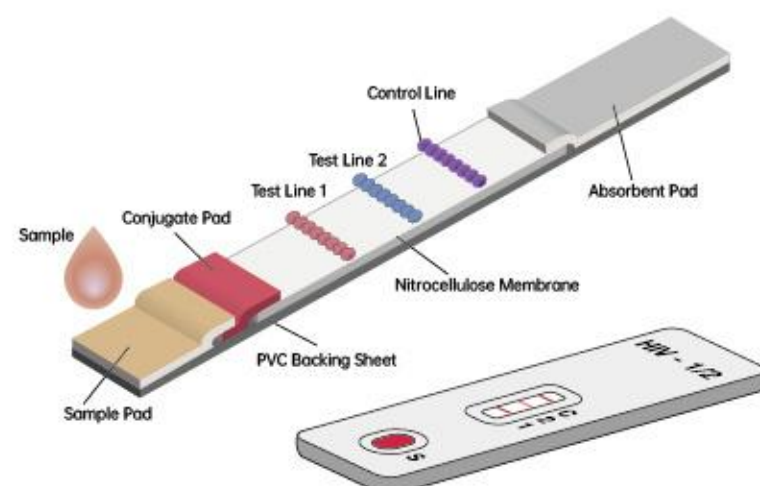

Advantages and Competitiveness



Client-Friendly · User-Friendly

Sensitivity · Specificity · Stability



Customizability · Flexibility

Cost-Effective (OEM · Private Labeling)



Comprehensive · Competitive

Quantitative · Semi-Quantitative · Qualitative



Whole Process Management System


R & D · QC · BD (Global)



Professional Service Team

Manufacture · Delivery

Production Process Flowchart of POCT Finished Products

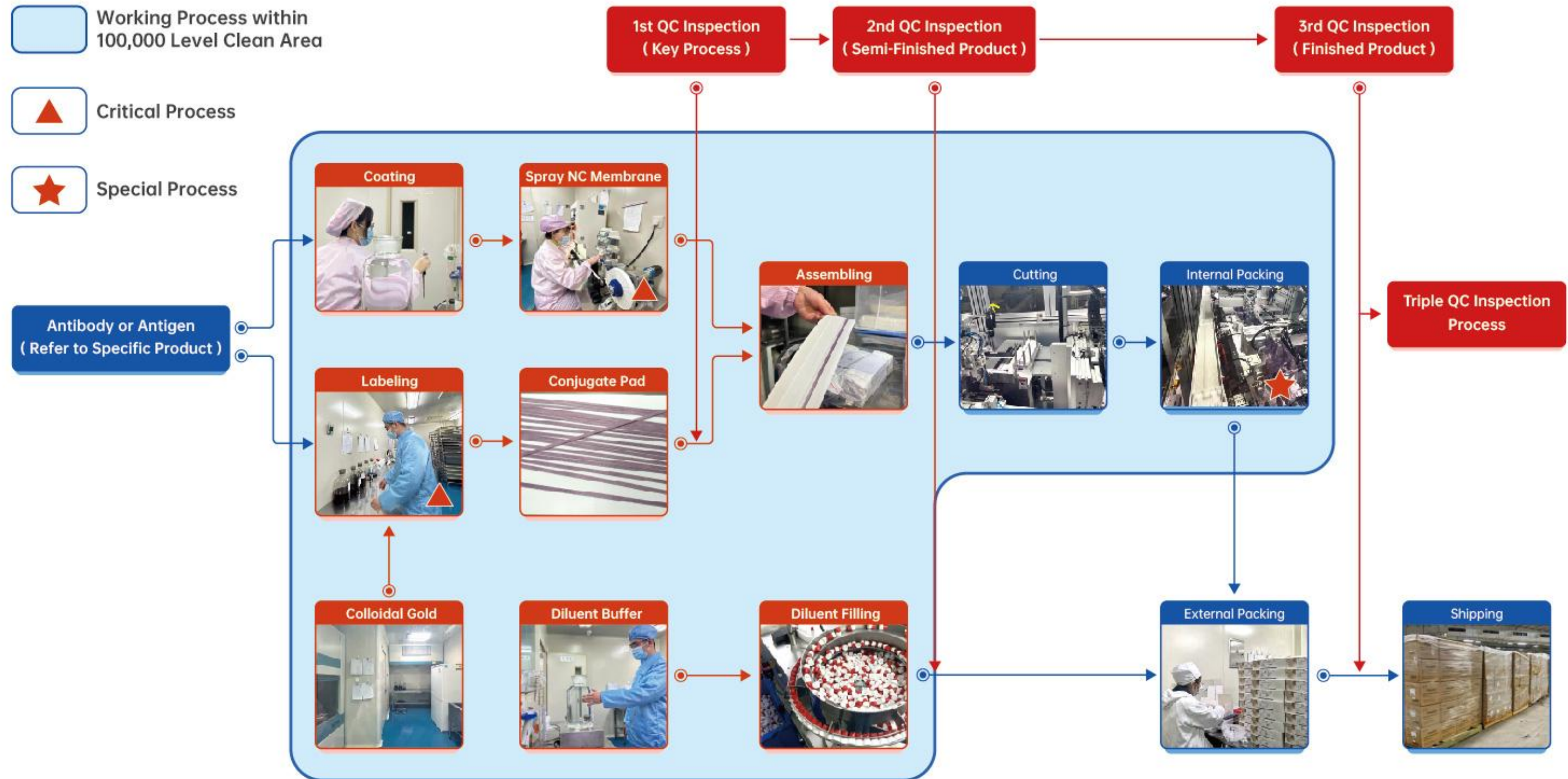

Remark : Uncut sheet (Semi - Finished Product) is available for shipping from NewScen after 2nd QC Inspection for further operation and processing from client.

Category of Uncut Sheet

INFECTIOUS DISEASE TEST KIT

	UNCUT SHEET	PRODUCT DESCRIPTION	SPECIMEN	SPECIFICATION
★	HIV 1+2	3 Lines : HIV (1+2) Antibody Rapid Test Kit	WHOLE BLOOD / SERUM / PLASMA	Cassette : 300 * 70 / 60 / 65 mm Strip / Dipstick : 300 * 80 / 70 mm
	HIV 1/2	2 Lines : HIV (1/2) Antibody Rapid Test Kit		
★	HCV	HCV Antibody Rapid Test Kit		
★	HBsAg	Hepatitis B Surface Antigen (HBsAg) Diagnostic Test Kit		
	HBsAb	Hepatitis B Surface Antibody (HBsAb) Diagnostic Test Kit		
	HEV	Hepatitis E Antibody (HEV IgM) Diagnostic Test Kit		
★	T P	Treponema Pallidum (TP) Antibody Rapid Test Kit		
	T B	M. Tuberculosis (TB IgG) Antibody Rapid Test Kit	SERUM / PLASMA (WHOLE BLOOD NOT APPLIED)	Cassette : 300 * 70 / 60 / 65 mm
	HAV	Hepatitis A Antibody (HAV IgM) Diagnostic Test Kit		
	HPAb	H. Pylori (HP) Antibody Rapid Test Kit	WHOLE BLOOD / SERUM / PLASMA	
★	HPAg	H. Pylori (HP) Antigen Rapid Test Kit	STOOL	
	ADENOVIRUS (Ag)	Adenovirus Antigen Rapid Test Kit		
	NORO (Ag)	Norovirus (NORO) Antigen Rapid Test Kit		
	ROTA (Ag)	Rotavirus (ROTA) Antigen Rapid Test Kit		
	TYPHOID (Ag)	Typhoid Antigen Rapid Test Kit		
	TYPHOID (IgG/IgM)	Typhoid IgG/IgM Antiobdy Rapid Test Kit	WHOLE BLOOD / SERUM	
	HSV (IgG)	Herpes Simplex Virus (HSV Type II) IgG Rapid Test Kit	SERUM / PLASMA (WHOLE BLOOD NOT APPLIED)	
	RSV (Ag)	Respiratory Syncytial Virus (RSV) Antigen Rapid Test Kit	NASOPHARYNGEAL SWAB	
	STREP A (Ag)	Group A Streptococcal Antigen (STREP A) Antigen Rapid Test Kit	SECRETION	
	CHLAMYDIA (Ag)	Chlamydia Antigen Rapid Test Kit		

TROPICAL (MOSQUITO TRANSMITTED) DISEASE TEST KIT

UNCUT SHEET	PRODUCT DESCRIPTION	SPECIMEN	SPECIFICATION
MAR (Pf/Pan)	Malaria Pf/Pan Antigen Rapid Test Kit	WHOLE BLOOD (SERUM / PLASMA NOT APPLIED)	Cassette : 300 * 70 / 60 / 65 mm
★ DENGUE (IgG/IgM)	Dengue IgG/IgM Rapid Test Kit	WHOLE BLOOD / SERUM / PLASMA	
★ DENGUE (NS1 Ag)	Dengue NS1 Antigen Rapid Test Kit		
ZIKA (IgG/IgM)	Zika Virus IgG/IgM Rapid Test Kit		
ZIKA (NS1 Ag)	Zika NS1 Antigen Rapid Test Kit		
CHIK (IgG/IgM)	Chikungunya IgG/IgM Rapid Test Kit		

DRUGS OF ABUSE (DOA) TEST KIT

UNCUT SHEET	PRODUCT DESCRIPTION	SPECIMEN	SPECIFICATION
5-IN-1	One Step (MOR / MET / KET / THC / MDMA) Multiple Drug Test Kit	URINE	Cassette : 300 * 70 / 60 / 65 mm Strip / Dipstick : 300 * 80 / 70 mm
★ MOR	Morphine (MOR) Gold Rapid Screen Test Kit		
★ MET	Methamphetamine (MET) Gold Rapid Screen Test Kit		
KET	Ketamine (KET) Rapid Test Kit		
THC	Tetrahydrocannabinol (THC) Rapid Test Kit		
MDMA	Methylenedioxymethamphetamine (MDMA) Rapid Test Kit		
AMP	Amphetamine (AMP) Rapid Test Kit		
BZD	Benzodiazepine (BZD) Rapid Test Kit		
COC	Cocaine (COC) Rapid Test Kit		
MTD	Methadone (MTD) Rapid Test Kit		

TUMOR MARKER TEST KIT

UNCUT SHEET	PRODUCT DESCRIPTION	SPECIMEN	SPECIFICATION
PSA	Prostate Specific Antigen (PSA) Diagnostic Test Kit	WHOLE BLOOD / SERUM / PLASMA	Cassette : 300 * 70 / 60 / 65 mm
AFP	Alpha Fetoprotein (AFP) Quantitative Rapid Test Kit		
CEA	Carcino-Embryonic Antigen Quantitative Rapid Test Kit		
★ FOB	One-Step Fecal Occult Blood (FOB) Diagnostic Test Kit	STOOL	(FOB Strip : 300 * 80 / 70 mm)
T F	Transferrin (TF) Rapid Test Kit		

RESPIRATORY DISEASE (COVID-19 / INFLUENZA) TEST KIT

UNCUT SHEET	PRODUCT DESCRIPTION	SPECIMEN	SPECIFICATION
COVID-19 (Ag)	COVID-19 Antigen Rapid Test Cassette	THROAT / NASAL SWAB	Cassette : 300 * 70 / 60 / 65 mm
COVID-19 Ag & Flu (A+B)	COVID-19 Ag. & Influenza Combo Rapid Test Kit		
Flu (A/B)	Influenza (A/B) Rapid Test Kit		
Flu (A+B)	Influenza (A+B) Combo Rapid Test Kit		
COVID-19 (IgG/IgM)	3 Lines : COVID-19 IgG/IgM Rapid Test Cassette	WHOLE BLOOD / SERUM / PLASMA	
COVID-19 (NAb)	COVID-19 Neutralizing Antibody Rapid Test Cassette		

HORMONE TEST KIT

UNCUT SHEET	PRODUCT DESCRIPTION	SPECIMEN	SPECIFICATION
HCG	Human Chorionic Gonadotropin (hCG) Tests (Pregnancy tests)	URINE	Midstream : 300 * 70 / 60 / 65 mm
L H	Luteinizing Hormone (LH) Tests (Ovulation tests)		Cassette : 300 * 70 / 60 / 65 mm
FSH	Follicle Stimulating Hormone (FSH) Tests (Menopause Tests)		
TSH	Human Thyroid Stimulating Hormone (TSH) Quantitative Rapid Test Kit	WHOLE BLOOD / SERUM / PLASMA	Cassette : 300 * 70 / 60 / 65 mm

IRON DEFICIENCY TEST KIT

UNCUT SHEET	PRODUCT DESCRIPTION	SPECIMEN	SPECIFICATION
★ S F	Serum Ferritin (SF) Rapid Test Kit	WHOLE BLOOD / SERUM / PLASMA	Cassette : 300 * 70 / 60 / 65 mm

CARDIAC MARKER TEST KIT

UNCUT SHEET	PRODUCT DESCRIPTION	SPECIMEN	SPECIFICATION
D-Dimer	D-Dimer Quantitative Rapid Test Kit	WHOLE BLOOD / PLASMA	Cassette : 300 * 70 / 60 / 65 mm
cTnI	Cardiac Troponin-I (cTnI) Quantitative Rapid Test Kit	WHOLE BLOOD / SERUM / PLASMA	

INFLAMMATORY MARKER TEST KIT

UNCUT SHEET	PRODUCT DESCRIPTION	SPECIMEN	SPECIFICATION
CRP	C-Reactive Protein (CRP) Quantitative Rapid Test Kit	WHOLE BLOOD / SERUM / PLASMA	Cassette : 300 * 70 / 60 / 65 mm
PCT	Human Procalcitonin (PCT) Quantitative Rapid Test Kit		


RENAL INJURY TEST KIT

UNCUT SHEET	PRODUCT DESCRIPTION	SPECIMEN	SPECIFICATION
mAlb	Urine Microalbumin (mAlb) Quantitative Rapid Test Kit	URINE	Cassette : 300 * 70 / 60 / 65 mm



Accessory

- Disposable Plastic Dropper (In Bag)
- Sample Diluent Buffer (In Bottle)
- Cassette / Midstream (Empty Housing)
- Sterile Safety Lancet (CE Mark)
- Throat / Nasal Swab (CE Mark)
- Sample Collection Tube (Secretion)
- Alcohol Pad (CE Mark)
- Aluminum Foil Pouch (Customized)
- Biohazardous Waste Bag (Customized)
- Sample Collection Tube (Stool)
- Desiccant (CE Mark)
- Pipette
- Centrifuge

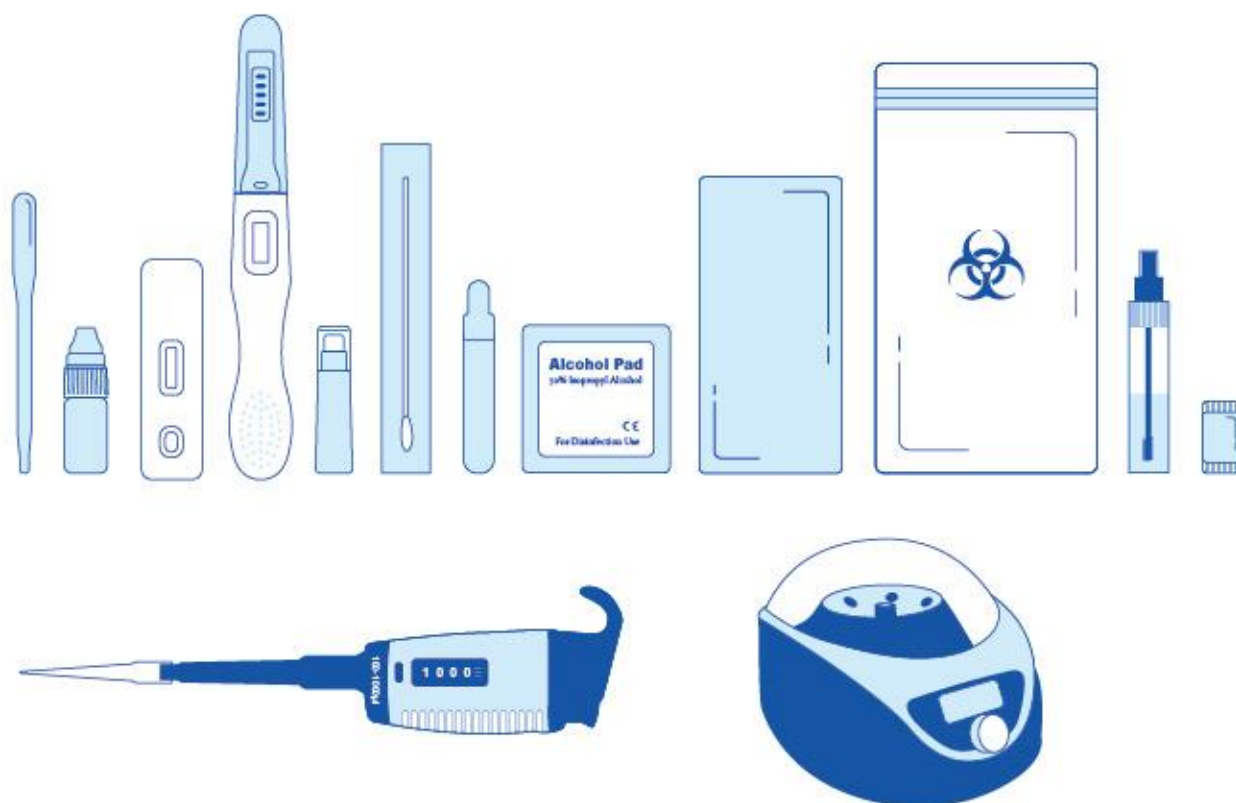

★ Remark : Accessories are optional and flexible for OEM requirements.
Equipments are optional for clients and / or end-users.

FAQs

? Basic Info (NewScen · Manufacturer)

About us :

- MOQ : 200 Uncut Sheets
- Method : Lateral Flow Assay (Immunochromatography Assay)
- Capacity : \cong 10,000 Uncut Sheets (Regular Daily Output / Production Line)
- Category : 10 Categories (57 Kinds of Uncut Sheets in 10 Categories)

? Trial Order (OEM · Private Labeling)

About you :

- Standard : QC Inspection Standard (Preferably in a Check Form, or the descriptions in details including requirements for the Performance / Appearance, etc.)
- Design : Design and Dimensions Drawing with precise sheet component position in details
- Format : Cassette / Midstream / Strip (Preferably in Picture with Dimension)

? Sample (Regular Format · Evaluation)

For uncut sheet sample from NewScen, please contact us by website / e-mail / social media :

- Website : www.newscentestkit.com
- E-mail : export@newscen.com
- Social Media : Facebook / Youtube / LinkedIn / Instagram / Twitter

